

**Date:** 02/03/2016

**Instruments:** *Convergys*<sup>®</sup> X5

**Subject:** Pressure sensor board part number revision

**Effective Date:** Immediate

Dear Customer,


We would like to inform you of an important new spare part number revision modification for the pressure sensor board on the *Convergys*<sup>®</sup> X5 instrument.

The new part number for the pressure sensor is now -A5519\_V3. This replaces P/N: A5519\_V2.

This change is necessary as we did not make it clear in the bulletin on 30<sup>th</sup> November 2015 (TPB NO-30112015) regarding the pressure sensor, that P/N A5519\_V2 was in fact the part number for a complete pressure sensor kit (i.e. consisting of a pressure sensor board, a puffer tank and tubing).

This means that it is now possible to purchase just the pressure sensor board rather than the whole kit. The pressure sensor board (A5519\_V3) must be used with puffer tank version (A510\_V2) and must be fitted with tubing of 14cm which is 2 cm longer as previously described.

### Pressure board Images

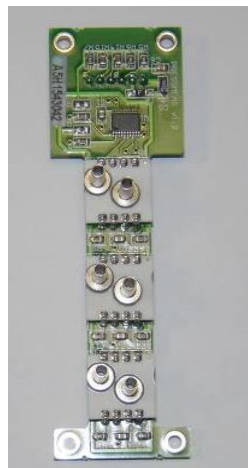

A5519v3

---

## Convergent Technologies GmbH & Co. KG

### **Registered Office**

Germany  
Ringstrasse 14, 35091 Cölbe, Germany  
Tel: +49 6421 8869948 Fax: +49 6421 8869958  
**Homepage:** [www.convergent-technologies.de](http://www.convergent-technologies.de)  
**e-mail:** [info@convergent-technologies.de](mailto:info@convergent-technologies.de)

**Ust-Id:** DE 241291597  
Steuernummer: 031 310 30121  
Amtsgericht Marburg HRA 4125  
**Geschäftsführer:** Umesh Bhalla

**Bank:** Volksbank Main-Tauber eG, Tauberbischofsheim,

### **EURO Transfers**

Account No. 67131401  
IBAN DE66 6739 0000 0067 1314 01  
BIC GENODE61WTH  
SWIFT GENODEFF

### **USD Transfers**

Account No. 567131408  
IBAN DE97 6739 0000 0567 1314 08  
BIC GENODE61WTH  
SWIFT GENODEFF

REV006

For details of your specific pricing for this P/N –A5519\_V3 please consult your business manager.

We apologise for this error and any inconvenience caused but hope that this is now clear.

Best regards,

**Convergent Technologies GmbH & Co. KG**  
**Ringstrasse 14, 35091 Coelbe, Germany**  
**Tel: +49 6421 8869948**  
**Fax: +49 6421 8869958**

---

**Convergent Technologies GmbH & Co. KG**

**Registered Office**

Germany  
Ringstrasse 14, 35091 Cölbe, Germany  
Tel: +49 6421 8869948 Fax: +49 6421 8869958  
**Homepage:** [www.convergent-technologies.de](http://www.convergent-technologies.de)  
**e-mail:** [info@convergent-technologies.de](mailto:info@convergent-technologies.de)

**Ust-Id:** DE 241291597  
Steuernummer: 031 310 30121  
Amtsgericht Marburg HRA 4125  
**Geschäftsführer:** Umesh Bhalla

**Bank:** Volksbank Main-Tauber eG, Tauberbischofsheim,

**EURO Transfers**

Account No. 67131401  
IBAN DE66 6739 0000 0067 1314 01  
BIC GENODE61WTH  
SWIFT GENODEFF

**USD Transfers**

Account No. 567131408  
IBAN DE97 6739 0000 0567 1314 08  
BIC GENODE61WTH  
SWIFT GENODEFF

REV006