



Convergys® SARS-CoV-2 Rapid Antigen Test

for Acute COVID-19 Detection


The Convergys® COVID-19 Rapid Antigen Test is a rapid chromatographic immunoassay for the qualitative detection of the Nucleocapsid protein of SARS-CoV-2, the pathogen causing the COVID-19 disease, from nasopharyngeal patient swabs directly, or after swabs have been added to viral transport media.

It is intended to aid in the rapid diagnosis of acute COVID-19 infections. The test is fast, cost-efficient and simple, which makes it an ideal addition for direct and acute COVID-19 testing when real-time PCR is too slow, costly and laborious.

For professional use only!

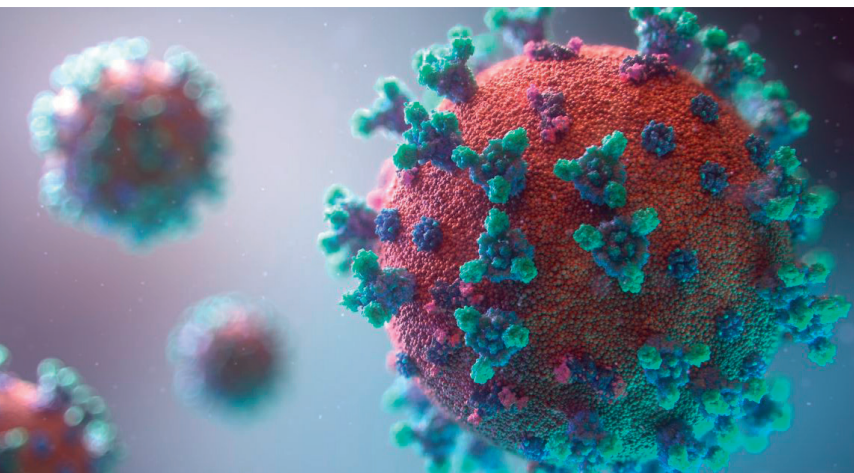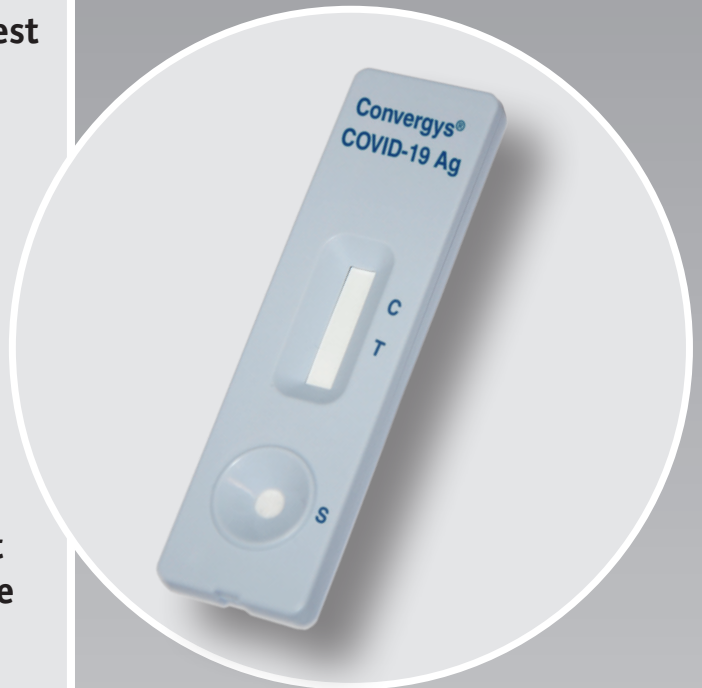

- Rapid testing for SARS-CoV-antigen within 15 minutes
- Simple and time-saving
- Includes all required materials, no further equipment necessary

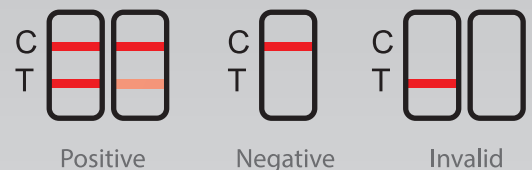

Specifications

Time to result:	15 minutes
Storage:	2-30°C
Shelf life:	24 months
Specimen type:	Nasopharyngeal swab
Relative Sensitivity:	96.7%
Relative Specificity:	99.2%
Accuracy:	98.7%

Contents

25 Test cassettes
25 Sterile swabs
25 Extraction tubes with dropper tips
Buffer solution for 25 tests
Tube rack
Package insert

Ordering information

REF	Article Name
1100-5330	Convergys® COVID-19 Rapid Antigen Test Kit